



FLIR A6750sc SLS

Longwave Infrared Thermal Camera

The FLIR A6750sc SLS incorporates a cooled Strained Layer Superlattice (SLS) detector that operates in the 7.5 to 9.5 micron waveband producing crisp LWIR thermal imagery at 640x512 pixel resolution.

FAST INTEGRATION TIMES

Working in snapshot mode, the FLIR A6750sc SLS is able to capture all pixels from a scene simultaneously in under 190 μ s for room temperature scenes. This is particularly important when monitoring fast moving objects where an uncooled thermal imaging camera would suffer from image blur. The camera supports image frame rates up to 4.1k frames per second when operating in windowing mode.

STANDARD VIDEO INTERFACES

Using a standard GigE Vision[®] interface to transmit full dynamic range digital video, and GenICam for camera control, the FLIR A6750sc SLS is a true “plug and play” thermal imaging camera. Additional interfaces include a BNC analog video output. The Gigabit Ethernet and analog video are simultaneously active yet independently controlled allowing greater flexibility for recording and display purposes.

CUSTOM COLD FILTERS AVAILABLE

Custom cold filtering options for specific spectral detection and measurement are available.

SOFTWARE

FLIR A6750sc SLS camera works seamlessly with FLIR ResearchIR Max software enabling intuitive viewing, recording and advanced processing of the thermal data provided by the camera. A Software Developers Kit (SDK) is optionally available.

COMPATIBLE WITH 3RD PARTY SOFTWARE

Control the A6750sc SLS and capture data directly into MathWorks[®] MATLAB software for custom image analysis and enhancement.

KEY FEATURES

- EXCELLENT LWIR IMAGE QUALITY: 640 X 512 PIXELS
- HIGH SPEED IMAGE ACQUISITION:
UP TO 4.1kHz IN WINDOWING MODE
- SYNCHRONIZATION WITH OTHER INSTRUMENTS AND EVENTS
- WIDE TEMPERATURE RANGES UP TO +2000°C
- MATLAB COMPATIBILITY



Circuit board




Motorcycle disc brake system



Imaging Specifications

Detector		A6750sc SLS	
Detector Type	Strained Layer Superlattice (SLS)		
Spectral Range	7.5 - 9.5 μm		
Resolution	640 x 512		
Detector Pitch	15 μm		
NETD	<30mK		
Well Capacity	7.2 M electrons		
Operability	>98%		
Sensor Cooling	Closed Cycle Rotary		
Electronics / Imaging			
Readout	Snapshot (FLIR 4 Channel)		
Readout Modes	Asynchronous Integrate While Read; Asynchronous Integrate Then Read		
Synchronization Modes	Sync In, Sync Out, Trigger In		
Integration Time	480 ns to 687 sec		
Frame Rate (Full Window)	Programmable 0.0015Hz to 125Hz		
Subwindow Mode	User Defined Size, Centered in Image		
Max Frame Rate (@ Min Window)	4,175Hz (16 x 4)		
Dynamic Range	14-bit		
Digital Data Streaming	Gigabit Ethernet		
Analog Video	NTSC, PAL		
Command & Control	Gigabit Ethernet and RS-232		
Measurement			
Standard Temperature Range	-20°C to 650°C (-4°F to 1,202°F)		
Optional Temperature Range	Up to 1,500°C (2,732°F) Up to 2,000°C (3,632°F)		
Accuracy	$\pm 2^\circ\text{C}$ or $\pm 2\%$ of reading		
Optics			
Camera f/#	2.5 or 4.0		
Available Lenses	13mm, 25mm, 50mm, 100mm, 200mm		
Focus	Manual		
Filtering	Behind the Lens, Custom Cold Filtering		
Image Presentation			
Analog Palettes	Grayscale + Color		
AGC	Manual, Linear, Plateau Equalization, DDE		
Zoom	Video Zoom is Auto Selected: Full Res = 1x, 1/4 Res = 2x		
General			
Operating Temperature Range	-40°C to 50°C (-40°F to 122°F)		
Storage Temperature Range	-55°C to 80°C (-67°F to 176°F)		
Altitude	0 to 10,000 Feet Operational; 0 to 70,000 Feet Non-Operational		
Shock / Vibration	40 g , 11 msec 1/2 sine pulse / 4.3 g RMS Random Vibration, All 3 Axis		
Power	24 VDC (< 50 W steady state)		
Weight w/o Lens	5 lbs		
Size (L x W x H) w/o Lens	7.7" x 4.0" x 4.0"		
Mounting	2 x 1/4"-20, 1 x 3/8"-16, 4 x 10/24		

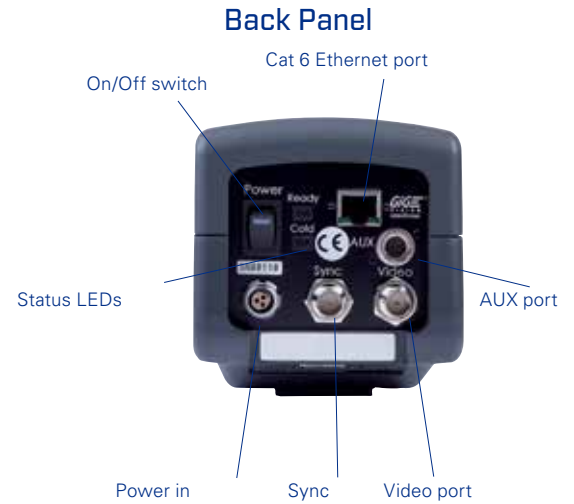

PORTLAND
Corporate Headquarters
FLIR Systems, Inc.
27700 SW Parkway Ave.
Wilsonville, OR 97070
USA
PH: +1 866.477.3687

BELGIUM
FLIR Systems
Luxemburgstraat 2
2321 Meer
Belgium
PH: +32 (0) 3665 5100

www.flir.com
NASDAQ: FLIR

NASHUA
FLIR Systems, Inc.
9 Townsend West
Nashua, NH 06063
USA
PH: +1 603.324.7611

UK
FLIR Systems UK
2 Kings Hill Avenue
Kings Hill
West Malling - Kent
ME19 4AQ
United Kingdom
PH: +44 (0) 1732 220 011

Equipment described herein may require US Government authorization for export purposes. Diversion contrary to US law is prohibited. Imagery for illustration purposes only. Specifications are subject to change without notice. ©2015 FLIR Systems, Inc. All rights reserved. (Updated 01/06/15)